

hawo

Seal**LINK** | Smart Networking System

# SIMPLE INTELLIGENT PACKAGING



## SMART NETWORKING SYSTEM

The **innovative** SealLINK networking system takes the packaging process for medical devices to a new, even more efficient level.


The direct interaction of a hawo sealing machine equipped with SealLINK software and our hm 8000 AS-V pouch production machine allows automatic and customised pouch production with just a single scan. Thanks to the intelligent SealLINK technology, pouches that are too short or too long are a thing of the past. The process is continuously optimised with each new instrument scanned by the data matrix reader.



SealLINK  
READY

The following models are compatible with SealLINK:


**hm 950 DC-V** NanoPak  
**hm 4000 DC-VI**  
**hm 5000 DC-VI** PrintPak

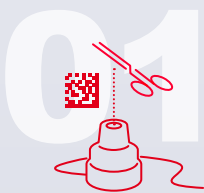


Scan the instrument code with the data matrix reader.



SealLINK automatically determines the appropriate pouch length and stores it.

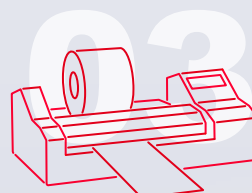

hm 8000 AS-V produces the right pouch directly with the next scan.

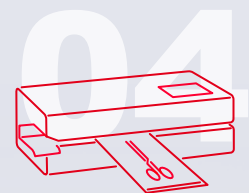

Insert the instrument and seal. Done!



Customised pouches up to max. 2,500 mm length.

SealLINK combines work efficiency, process optimisation and material savings in one single system.

**Scan one time, store data one time, get the right pouch size each and every time.**