

hawa

Seal**LINK** | Smart Networking System

SIMPLE INTELLIGENT PACKAGING



SMART NETWORKING SYSTEM


The **innovative** Seal**LINK** networking system takes the packaging process for medical devices to a new, even more efficient level.

The direct interaction of a hawo sealing machine equipped with SealLINK software and our hm 8000 AS-V pouch production machine allows automatic and customised pouch production with just a single scan. Thanks to the intelligent SealLINK technology, pouches that are too short or too long are a thing of the past. The process is continuously optimised with each new instrument scanned by the data matrix reader.



The following models are compatible with SealLINK:


hd 650 DET-V EcoPak 06
hd 650 DCT-V EcoPak 07
hm 950 DC-V NanoPak
hm 4000 DC-VI
hm 5000 DC-VI PrintPak

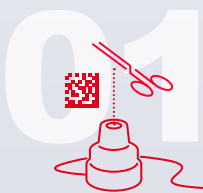


Scan the instrument code with the data matrix reader.



SealLINK automatically determines the appropriate pouch length and stores it.

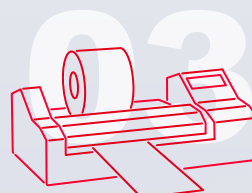

hm 8000 AS-V produces the right pouch directly with the next scan.

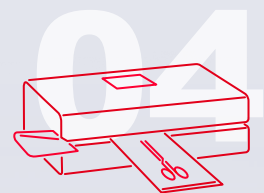

Insert the instrument and seal. Done!



Customised pouches up to max. 2,500 mm length.

Seal**LINK** combines work efficiency, process optimisation and material savings in one single system.

Scan one time, store data one time, get the right pouch size each and every time.